


# DUAL LED

Renew Your Skin Smart  
Light Therapy

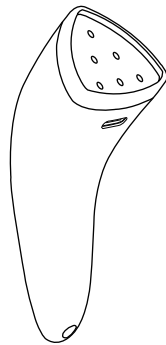


User Manual

## Description of the device

**PRECISE TARGETING LIGHT**  
Focuses on individual area, not disturbing the surrounding healthy skin

**POWER OF LIGHT THERAPY** The LED acne fitted with 6 high-intensity LED bulbs that provide professional skincare in the comfort of your own home Energy in the form of light waves are UV-free and non-thermal, absolutely safe on skin.

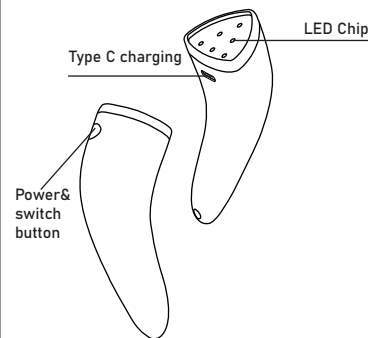

## How to use

1. Cleanse and dry your face before LED photon treatment.
2. Press "Start/Stop" button to turn on the device.
3. Hold the device in contact with the skin, slowly for 3 minutes each treatment
4. Use red light to treat wrinkles and fine lines. Slowly move on the treatment areas.
5. Use blue light to treat eczema and acne. Slowly move around the areas to aid in the natural healing
6. Apply skin care to moisturize your skin.

## Mode and display

Long press the power button for 2 seconds to turn on/off the device. Two gears can be selected. Press the power button to switch the gears.

--- First gear: Red Light

--- Second gear: Blue Light  
Low Battery: red light flashing  
Charging: red light flashing (breathing light)  
Fully charged: Solid blue light

## Attention

- 1, Please keep it out of the reach of infants.
- 2, Please do not use following symptoms: heart disease patients. Severe patients, acute patients, infectious patients, allergic physical patients, dermatitis, etc.
- 3, Please do not use during pregnancy and lactation.
- 4, Please do not use in the following areas: wound, plastic surgery, eyeball.
- 5, Keep away from fire and high temperature.
- 6, Please stop using the device when it is abnormal.

## Warning

- \*Please do not use other charging cables to avoid damaging the device.
- \*Please do not use disassemble or repair this product by yourself to avoid damage.
- \*Please avoid falling and hitting the product to avoid damage.
- \*Do not rinse the product directly into water.

## Product Parameters

Product Name	Renew Your Skin Smart Light Therapy
Material	ABS+LED chip
Working voltage	3.7V
Working current	0.06a
Working power	About 0.25w
Charging voltage	5V
Charging current	0.3a
Charging power	About 1.5w
Adapter requirements	5v/0.5a
Battery	250mAh
Red light	630NM±10NM
Blue light	450±10NM
Number of LEDs	6pcs
Automatic shutdown setting	3 minutes
Bare metal size	34.7*82*68.5mm
Bare metal weight	25G
Charging port	type C charging

## Cleaning and Maintenance

- 1, Wipe it off after use immediately.
- 2, This product can not be washed with water directly.
- 3, Do not use the sharp or corrosive or volatile diluent to clean the product.

## Environmental protection

- instructions
- 1, This product is in accordance with the regulations and of the management method for the restricted use of hazardous substances in the electrical and electronic products.
  - 2, Consumers will not have the leakage and precipitation of harmful substances during normal use, consumer can be used safely.

## Abandonment instructions

- 1, The battery must be removed from the appliance before abandoning the appliance.
- 2, When taking out the battery, the appliance must first cut off the power supply.
- 3, The battery should be disposed of safely.
- 4, Please pay attention to the heating surface on the device.